

FongAI Co., Ltd.

FongAI Assessment APP Operation Manual



The screenshot shows the login interface of the FongAI Assessment APP. At the top, there are two dropdown menus for 'Font Size' and 'English'. Below these is the company logo, a stylized green and blue figure, followed by the text '建豐健康科技' and 'FongAI Co., Ltd.' in Chinese and English respectively. Underneath is 'FongAI Assessment'. The login form includes an 'Account' field with 'fongai' entered, a 'Password' field with a toggle icon, and a 'Remember account' checkbox which is checked. A large green 'Login' button is at the bottom. The footer of the screen displays 'APP version / 1.1.0'.

Font Size English

 建豐健康科技
FongAI Co., Ltd.
FongAI Assessment

Account

fongai

Password

☒ Remember account

Login

APP version / 1.1.0



建豐健康科技
FongAI Co., Ltd.

TABLE OF CONTENTS

1	SHORT PHYSICAL PERFORMANCE BATTERY (SPPB) →
2	FONGAI ASSESSMENT DOWNLOAD PROCESS →
3	FONGAI ASSESSMENT OPERATION PROCESS →
4	DATA UPLOAD →
5	OTHER PRECAUTIONS →

SHORT PHYSICAL PERFORMANCE BATTERY (SPPB)

1. Balance Tests

The participant is asked to stand unsupported for 10 seconds with their feet in a certain position.

A. Side-by-side stand



1 point (Held for 10 seconds)



0 points (Less than 10 seconds)

B. Semi-tandem stand



1 point (Held for 10 seconds)



0 points (Less than 10 seconds)

C. Tandem stand



2 points (Held for 10 seconds)



1 point (Held for 3-9.99 seconds)



0 points (Held for <3 seconds)

Balance test score subtotal: _____

2. Gait Speed Test

The test measures gait speed. The participant walks to the line on the floor (3 or 4 meters) at a normal pace. The time is recorded.



For 4-Meter Walk:



4 points: < 4.82 sec



3 points: 4.82 to 6.20 sec



2 points: 6.21 to 8.70 sec



1 point : > 8.70 sec



0 points: unable to finish

3. Chair Stand Test

The time to stand up and sit down five times consecutively



4 points: <11.19 sec



3 points: 11.2–13.69 sec



2 points: 13.7–16.69 sec



1 point : 16.7–59.99 sec



0 points: > 60 sec

Total

/ 12 points


Download Process & Pre-operation Setup

Download: Please follow the process below to search or scan the QR code to download.

1. Scan the QR code or open the App Store, and search for "FongAI".

(Recommended device: requires iOS 15.2 or later)

SCAN
ME!

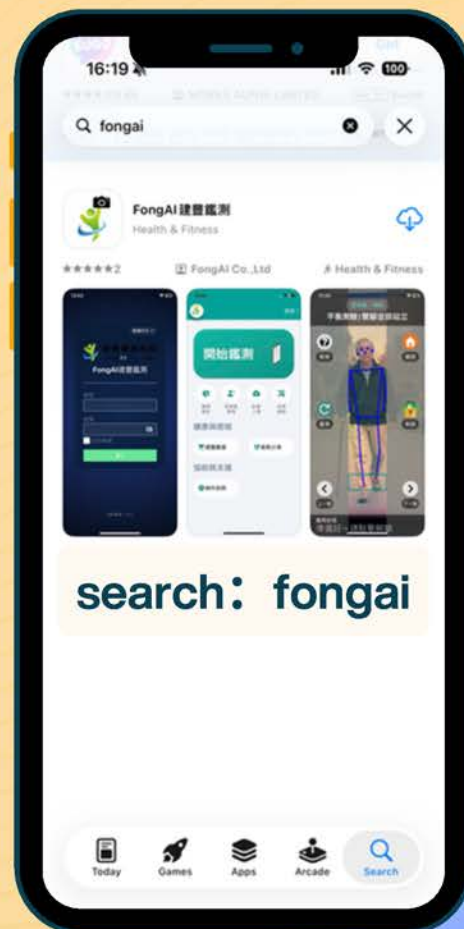

2. Download the "FongAI Assessment" APP



Assessment Instructions

The FongAI Dynamic Skeleton Fall-Prevention Assessment must be completed by two people working together:

one person as the participant and one person as the operator.

The assessment content is based on the Short Physical Performance Battery, also known as the SPPB. The participant will be assessed in the following order:

Chair Sit-to-Stand Test → Balance Test → Gait Speed Test



01.

Required Items



- One iPhone or iPad
- One lightweight tripod
- One stable chair



02.

Preparation



1. Place the iPhone or iPad vertically on the tripod. After fixing its position, do not move the tripod.
2. Place the chair in a stable position.
3. Make sure the phone and the chair are aligned in a straight line.
4. The camera lens should be approximately **95–97 cm** above the ground.
5. Adjust the device so that the **vertical angle is within 3 degrees** and the **horizontal level is within ± 0.03** .
6. Turn on the phone's sound. Voice prompts will be provided after the assessment begins.



03.

Camera Calibration

1. Open the FongAI Assessment App and log in.
2. For first-time use, follow the calibration prompt or tap “Camera Calibration.”
3. Enter the operator’s height, face the camera, and slowly step backward.
4. When the distance reaches 4 meters, you will hear **“Calibration Completed.”**

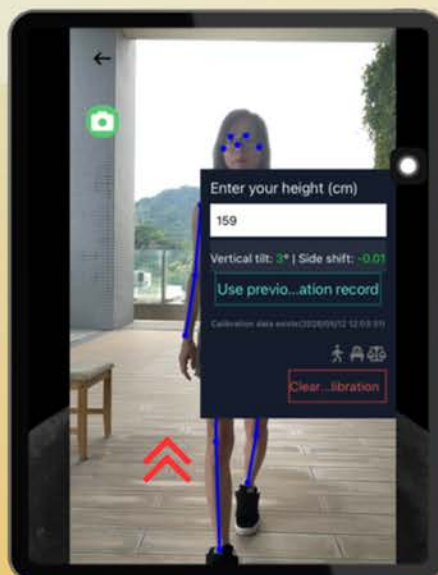
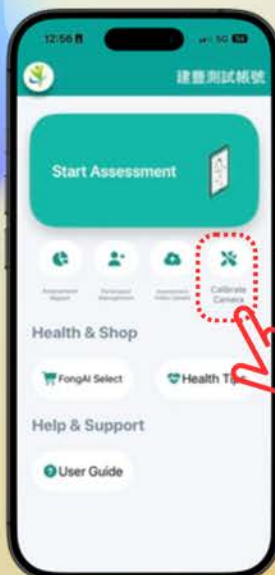


04.

Participant Management

1. On the App home page, tap **“Participant Management.”**
2. After entering the screen, select the folder, class, or unit where the participant list will be created.
3. Enter the participant’s name, gender, and date of birth.
4. After completing the information, tap Register.


***If you need to create a participant list in large quantities, please contact FongAI staff.**

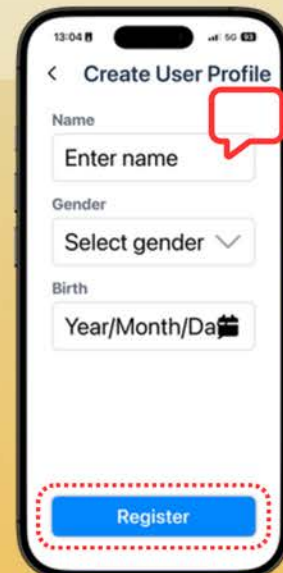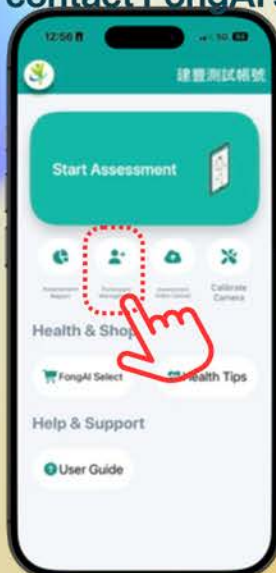

05.

Start Assessment



1. On the App home page, tap **"Start Assessment."**
2. After entering the screen, select the assessment site.
3. Swipe to select the participant who will be assessed, or use the magnifying glass icon in the upper-right corner to search the participant list.
4. Tap the arrow on the right to begin the assessment.

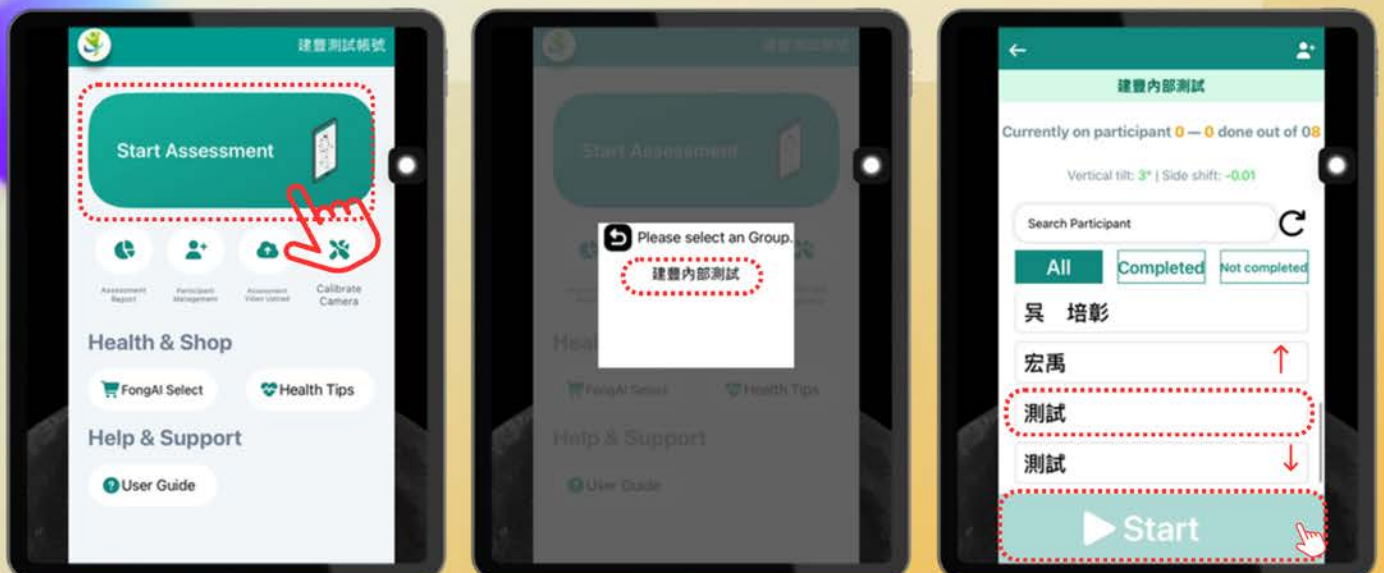

06.

Chair Sit-to-Stand Test


1. First, perform the sit-to-stand test. The participant should sit on a stable chair and place the ankle joint area within the green frame.
2. After the participant sits in the designated position, tap the Unlock button on the right side.
3. At this time, ask the participant to cross both arms over the chest and complete five consecutive sit-to-stand movements.

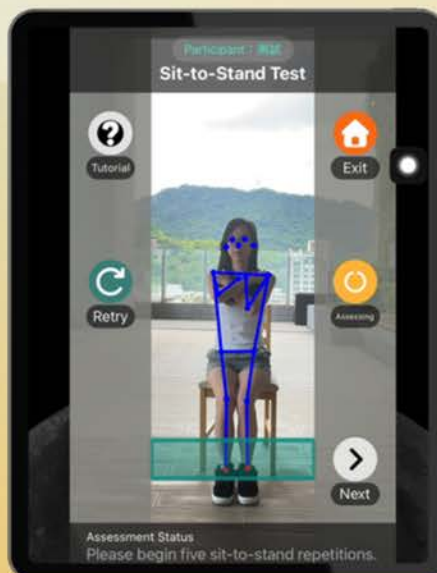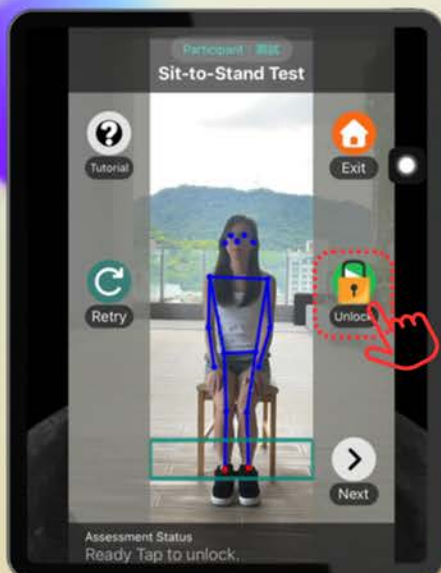

07.

Balance Test

1. Perform the balance test. Ask the participant to slowly stand up from the same position used in the sit-to-stand test, and move the chair away.
2. **Important: Do not allow the participant to sit down again.**
3. Align both ankle keypoints within the green frame.
4. The balance test includes:
 - Side-by-side | Semi-tandem | Tandem stand

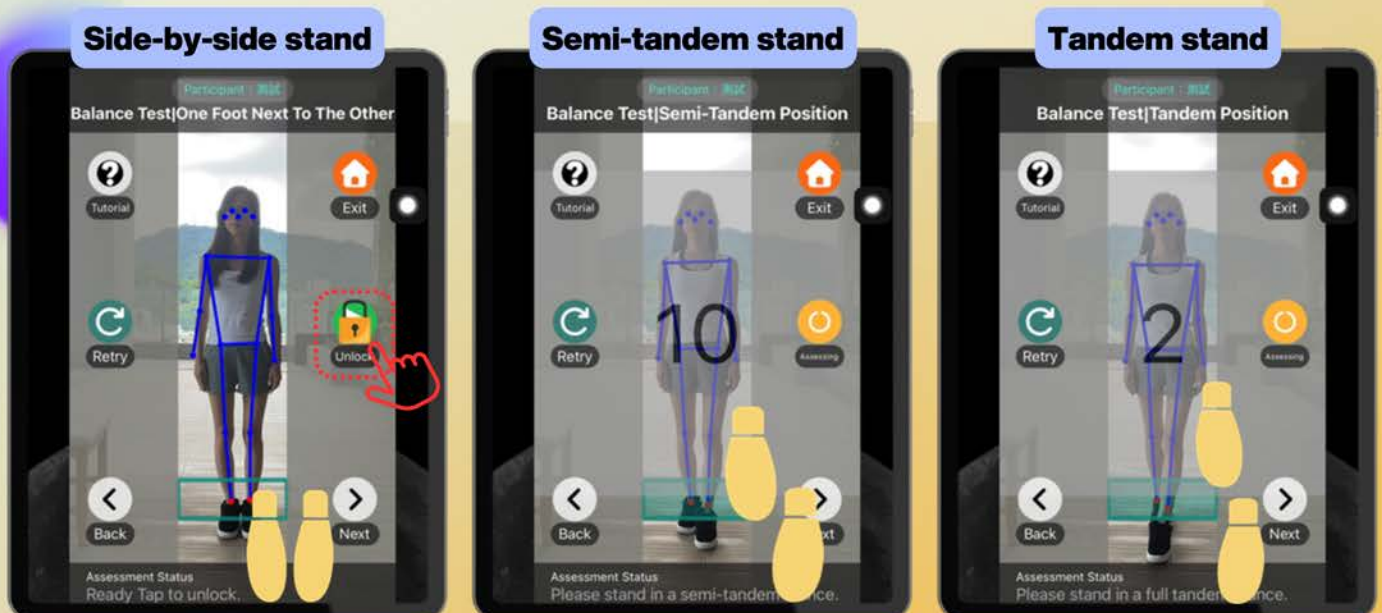

08.

Gait Speed Test

1. Perform the gait speed test. Ask the participant to stand in the orange area facing the device.
2. Tap Unlock and have the participant walk toward the device, then turn around after passing the green frame and return to the starting point.
3. The assessment is complete after finishing the walk test.

**Important: Ask the participant to walk at their usual pace.*




DATA UPLOAD

After completing the assessment, the video files temporarily stored on the phone or tablet must be uploaded in order to generate dashboard data.

The data upload process is as follows:

1. Return to the home page and tap **“Video Upload.”**
2. Select the date of the video to be uploaded.
3. Select the assessment site for the video to be uploaded, then tap Start Upload.

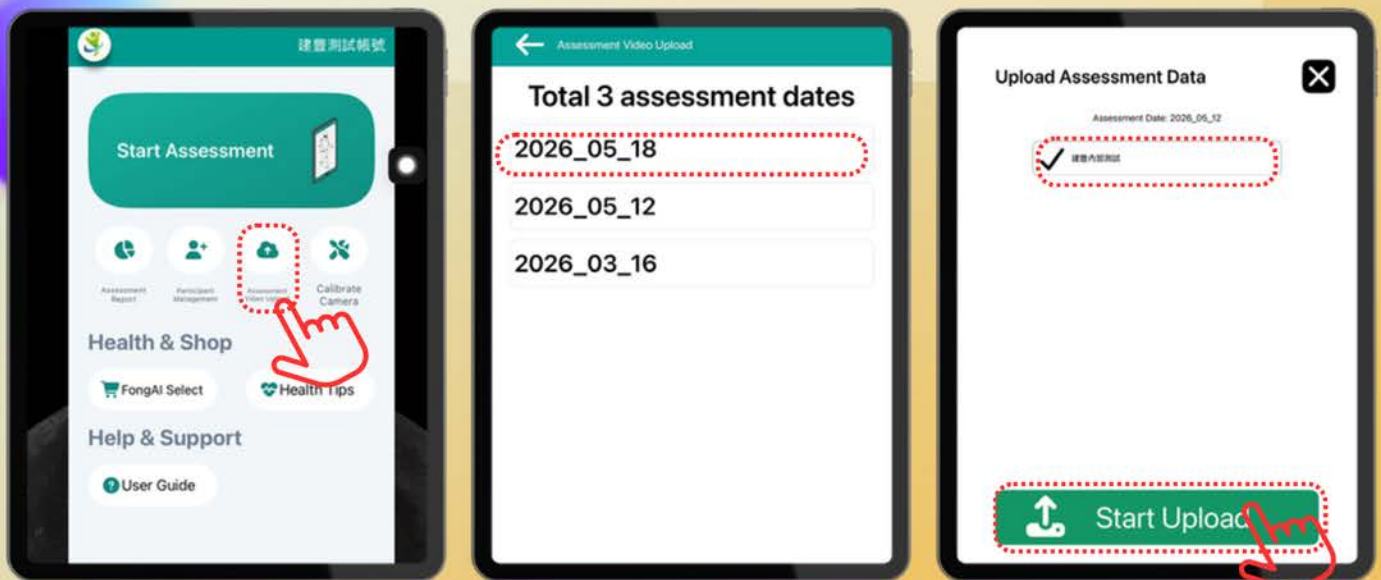

OTHER IMPORTANT NOTES

1.Button Descriptions

Re-test

If the older adult performs the movement incorrectly or is not ready, tap this button to restart the assessment.

Previous Step

If you feel the previous step was not completed properly, you may return to the previous assessment item.

Exit Assessment


Pause the assessment and return to the home page.

Unlock Button

Before each stage begins, first explain the movement to the older adult. After confirming that there are no issues, tap the Unlock button to begin.

Skip Button


If the older adult's physical condition prevents them from completing the stage, you may tap Skip.

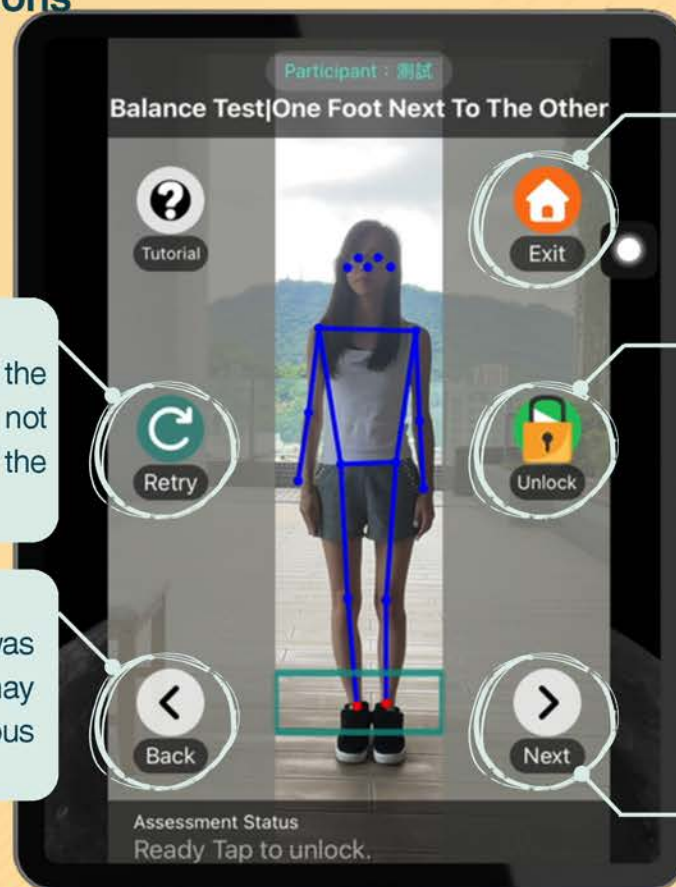

2.Manual Start

Manual Start

If the participant is wearing long pants, a skirt, or shoes of a similar color, the system may not be able to accurately detect the ankle keypoints. This may cause difficulty in automatically determining whether the movement is correct during the balance test. If the older adult's movement is correct, the operator may use Manual Start to allow the system to begin the countdown and conduct the assessment.

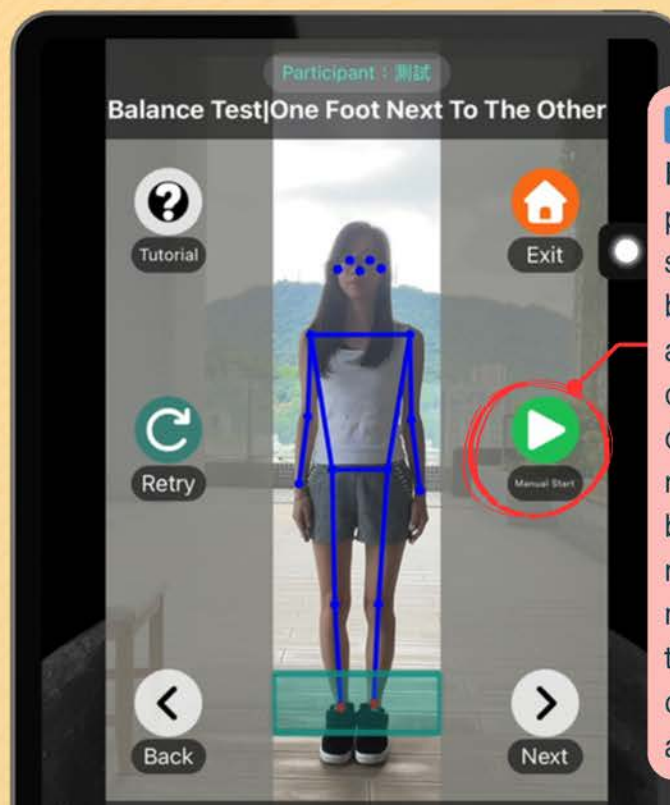

Important!! The Manual Start function can only be used during the balance test.

FongAI Co., Ltd.

Need Assistance?

If you encounter any system abnormalities, account issues, or process questions during operation, please feel free to contact us, and our specialists will assist you:

Support Email: service@fongai.co

Official Website: <https://fongai.co>

 **【 Safety Reminder 】** When performing the FongAI assessment, please strictly take **"the safety of the tested senior" as the highest guiding principle**. If the senior feels unwell or exhausted during the test, please stop the test immediately and do not force it.

Document Release Info.

Issuer: FongAI Co., Ltd.

Language: English (en-US)

Version: V1.1.6_EN

Update Date: 2026.05.28



建豐健康科技
FongAI Co., Ltd.